


AMIDES

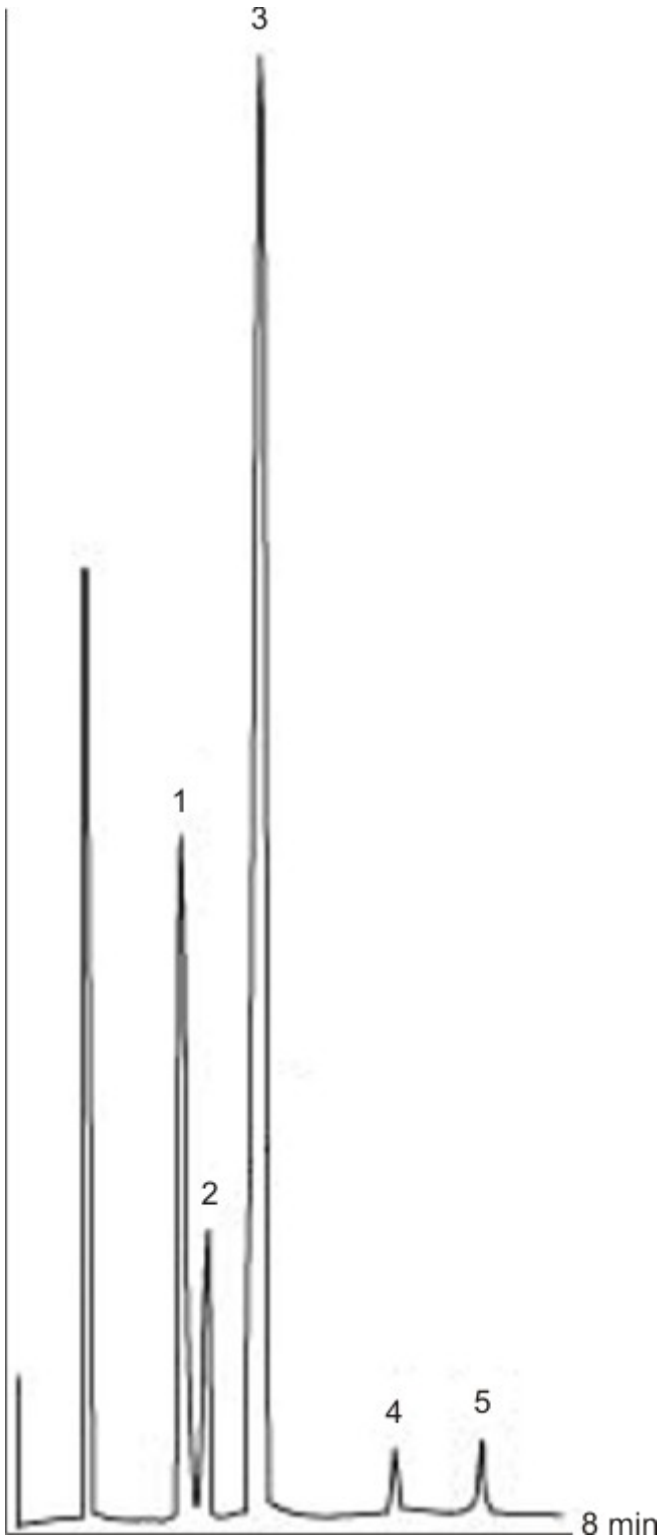

Column

Phase	MEGA-ACID
I.D.	0.53 mm
Film Thickness	1.0 µm
Length	10 m

Chromatographic Conditions

Injector	Split	
Injection Volume	1 µL	
Oven	T start	70°C
	Rate	10°C/min
	T end	200°C
Detector	FID	
Carrier Gas	Helium	4.0 mL/min

Peak Identification

1	3-Picoline
2	N,N-Dimethylformamide
3	N,N-Dimethylacetamide
4	N-Methylacetamide
5	Acetamide