


DIMETHYLANILINES

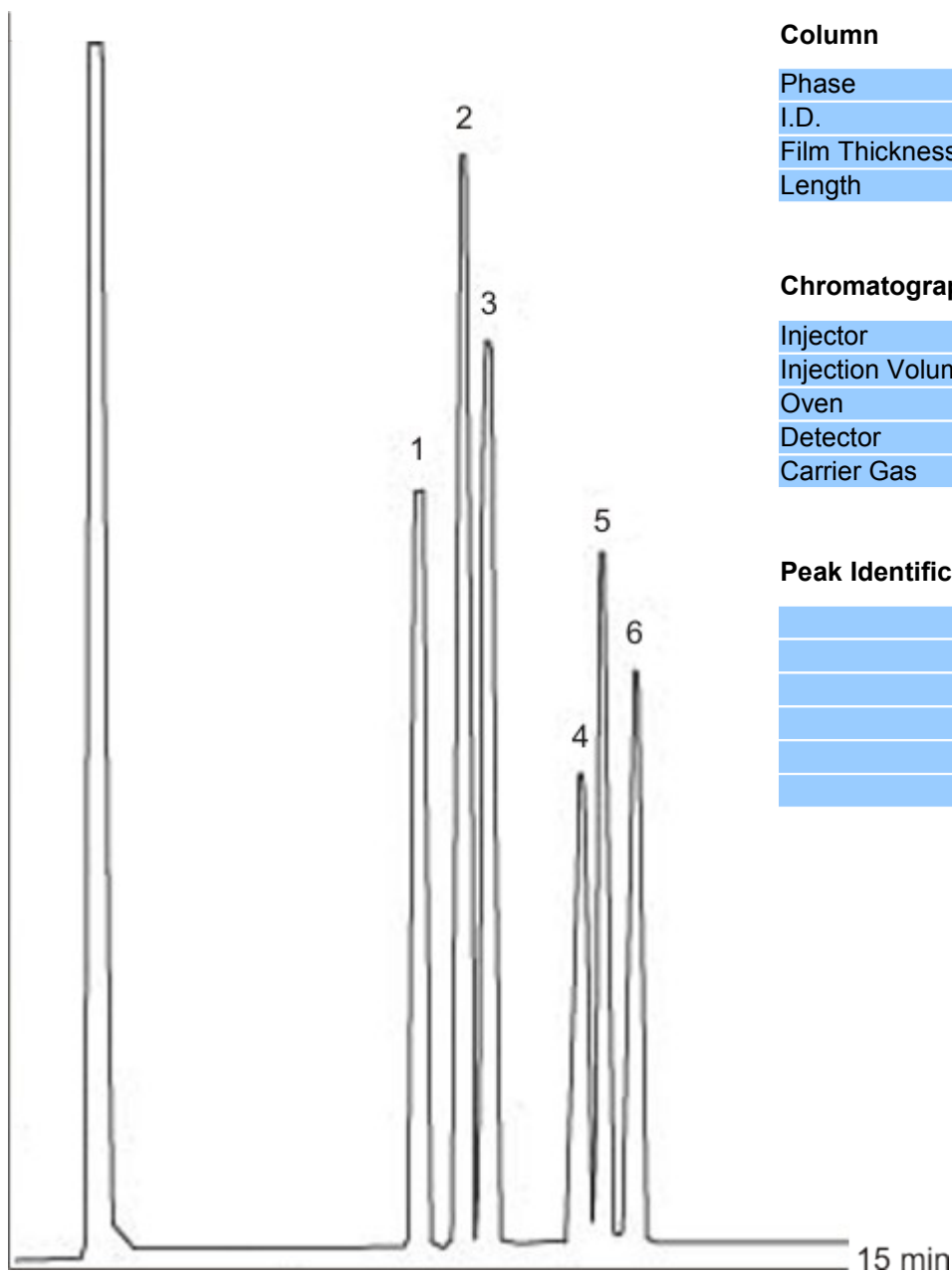

Column

Phase	MEGA-WAX
I.D.	0.32 mm
Film Thickness	0.25 µm
Length	30 m

Chromatographic Conditions

Injector	Split	
Injection Volume	1 µL	
Oven	Isot.	130°C
Detector	FID	
Carrier Gas	Nitrogen	1.8 mL/min

Peak Identification

1	2,6-Dimethylaniline
2	2,4-Dimethylaniline
3	2,5-Dimethylaniline
4	3,5-Dimethylaniline
5	2,3-Dimethylaniline
6	3,4-Dimethylaniline