

Bergamot Oil

Column: **MEGA-SE54 FAST - 0.10mm, 0.10 μ m, 5m**

Catalog Code: F-SE54-010-010-5

Conditions:


Injection: Split 230°C, 0.5 μ L, 1:250 Split Ratio.

Sample Dilution: 1% in Cyclohexane.

Detector: FID 250°C.

Oven Temp.: 50°C (0.1min), 30°C/min, 250°C.

Carrier Gas: Hydrogen, 0.5 mL/min.

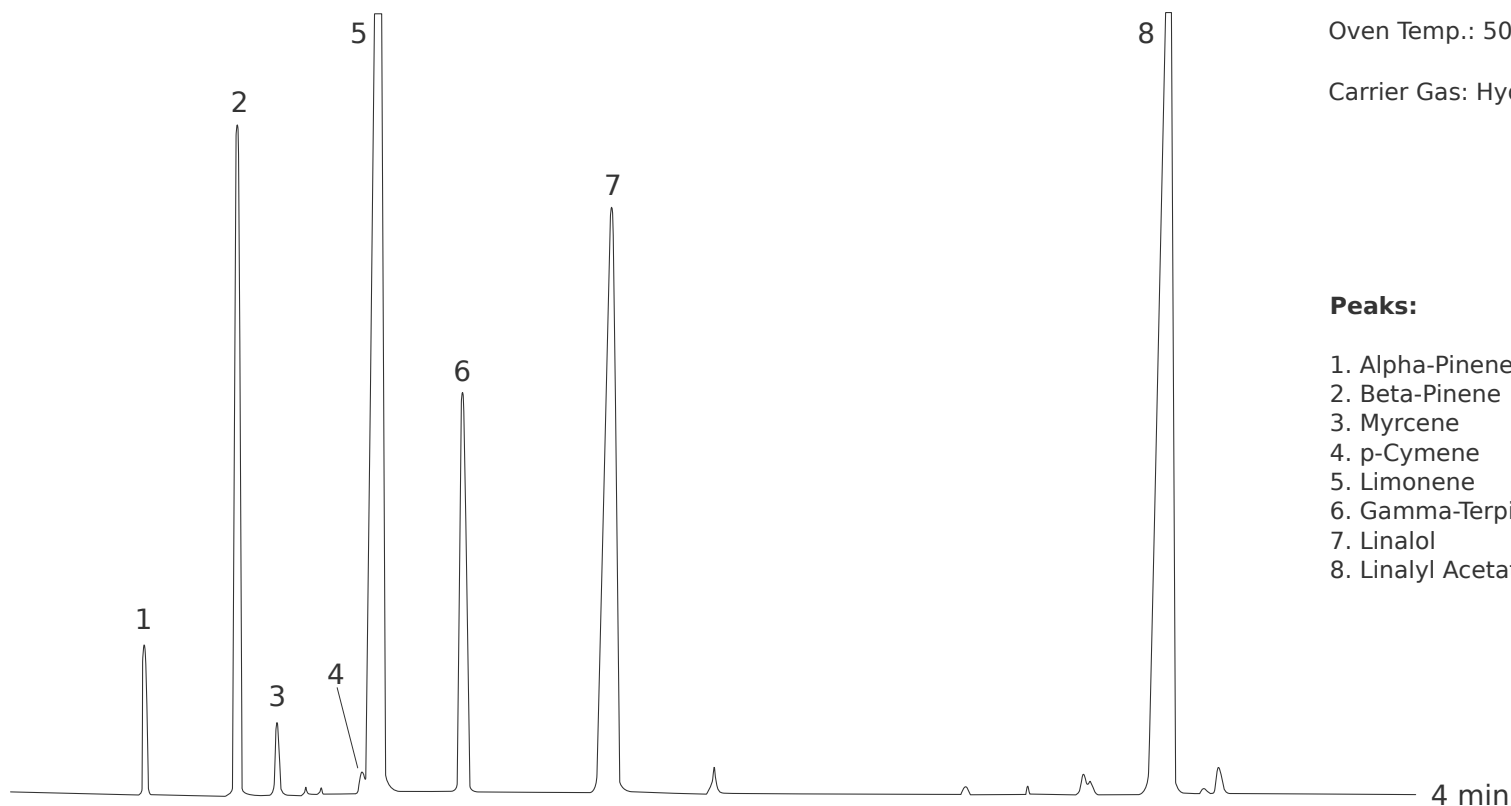

Peaks:

1. Alpha-Pinene
2. Beta-Pinene
3. Myrcene
4. p-Cymene
5. Limonene
6. Gamma-Terpinene
7. Linalol
8. Linalyl Acetate

Acknowledgement: Prof. C.Bicchi, C.Brunelli et al. - Università di Torino, Dipartimento di Scienza e Tecnologia del Farmaco, Via Pietro Giuria, 9 - Torino - Italy.

