


## FAMES C4 - C18:3 (Fatty Acid Methyl Esters)

Column: **MEGA-WAX - 0.32mm, 0.25 $\mu$ m, 25m**  
Catalog Code: S-WAX-032-025-25

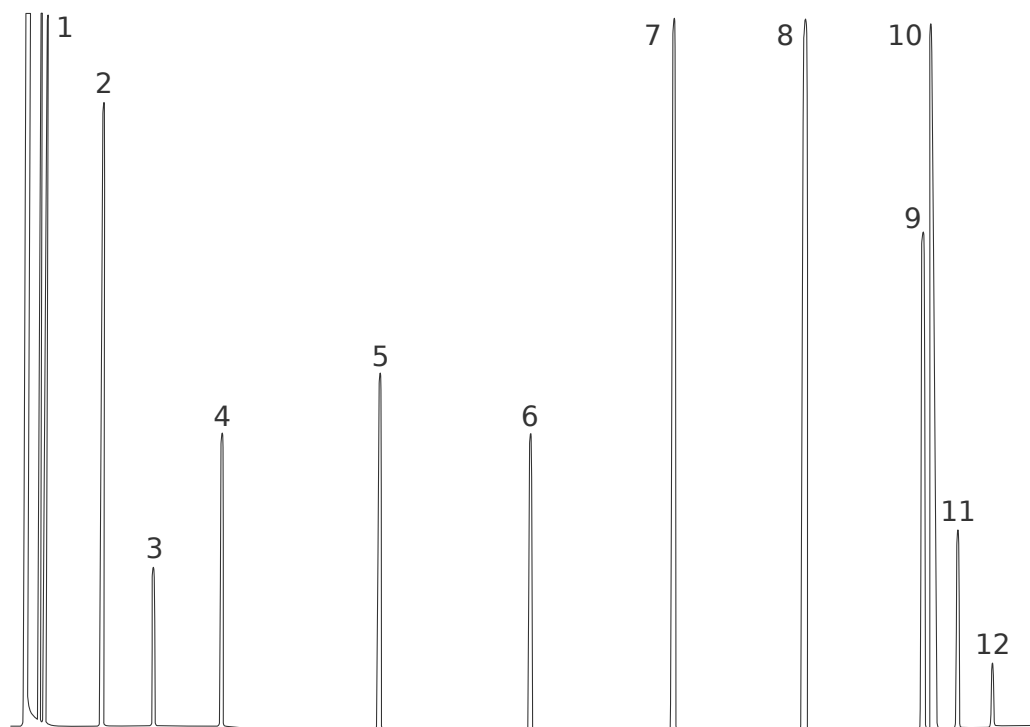

### Conditions:

Injection: Split 250°C, 1.0 $\mu$ L, 1:50 Split Ratio.

Detector: FID 250°C.

Oven Temp.: 50°C, 4°C/min, 210°C.

Carrier Gas: Hydrogen, 60 kPa.

### Peaks:

1. Butyric Acid ME (C4)
2. Caproic Acid ME (C6)
3. Heptanoic Acid ME (C7)
4. Caprylic Acid ME (C8)
5. Decanoic Acid ME (C10)
6. Lauric Acid ME (C12)
7. Myristic Acid ME (C14)
8. Palmitic Acid ME (C16)
9. Stearic Acid ME (C18:0)
10. Oleic Acid ME (C18:1)
11. Linoleic Acid ME (C18:2)
12. Linolenic Acid ME (C18:3)

