

Pesticides

Column: **MEGA-SE54 FAST - 0.10mm, 0.10µm, 5m**

Catalog Code: F-SE54-010-010-5

Conditions:

Injection: Split 230°C, 0.5µL, 1:200 Split Ratio.

Sample Dilution: 1% in Cyclohexane.


Detector: FID 250°C.

Oven Temp.: 50°C (0.1min), 40°C/min, 250°C.

Carrier Gas: Hydrogen, 0.5 mL/min.

Peaks:

1. Alpha-HCH
2. Gamma-HCH
3. Chlorothalonil
4. Heptachlor
5. Methyl Parathion
6. Ethyl Paraoxon
7. Malathion
8. Fenitrothion
9. Ethyl Parathion
10. /
11. Chlordane-*Trans*
12. Chlordane-*Cis* + Alpha-Endrin
13. Dieldrin
14. Beta-Endosulfan
15. o,p'-DDT
16. p,p'-DDT
17. Tetradifon

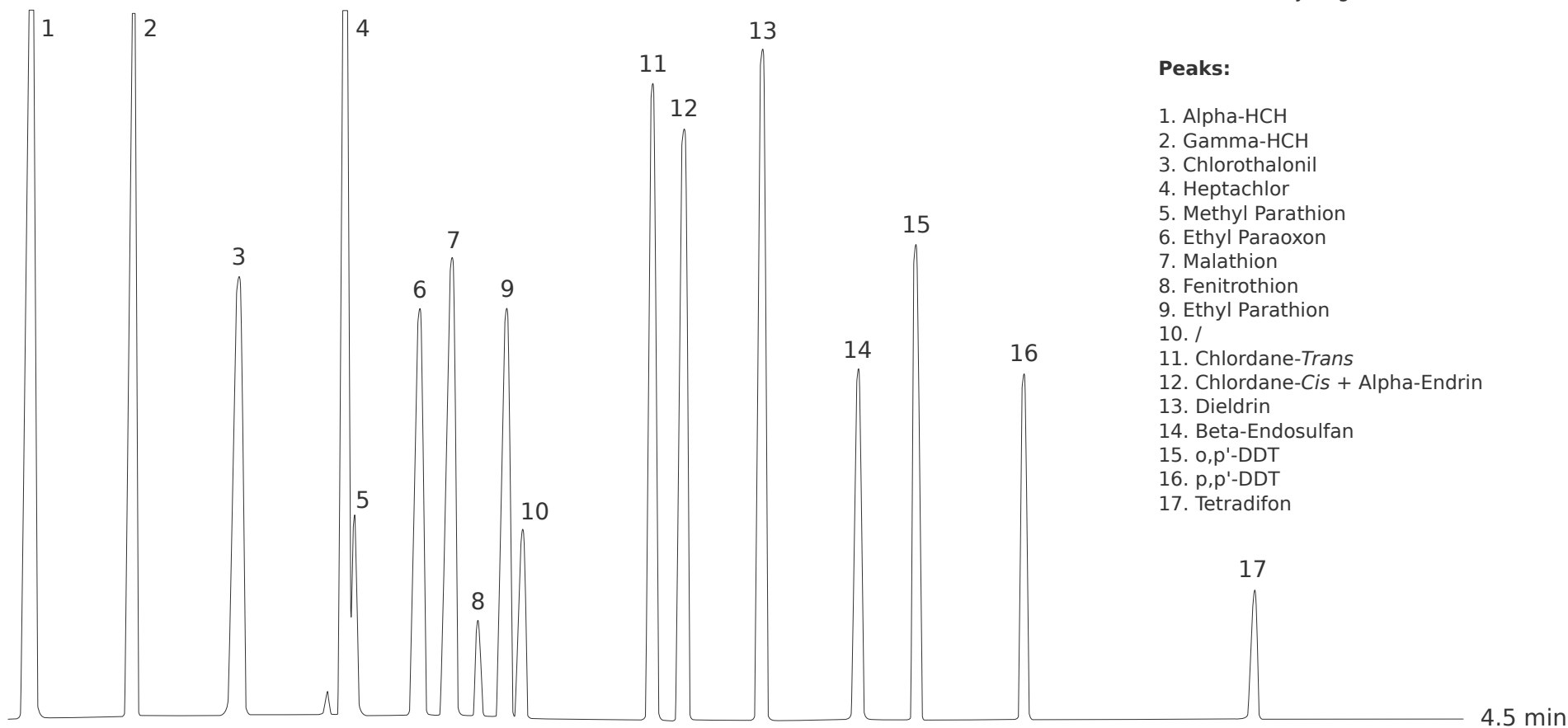

Acknowledgement: Prof. C.Bicchi, C.Brunelli et al. - Università di Torino, Dipartimento di Scienza e Tecnologia del Farmaco, Via Pietro Giuria, 9 - Torino - Italy.

